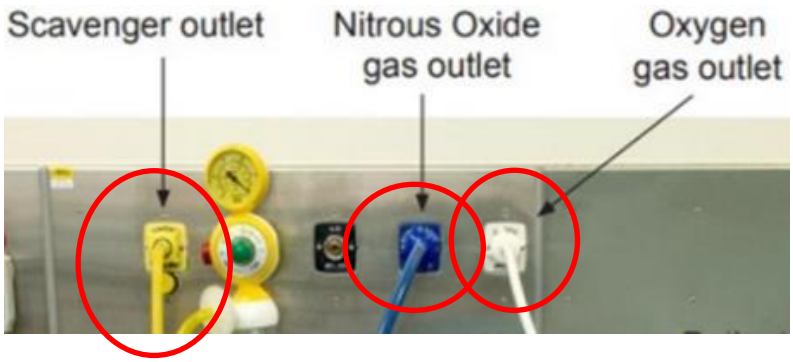
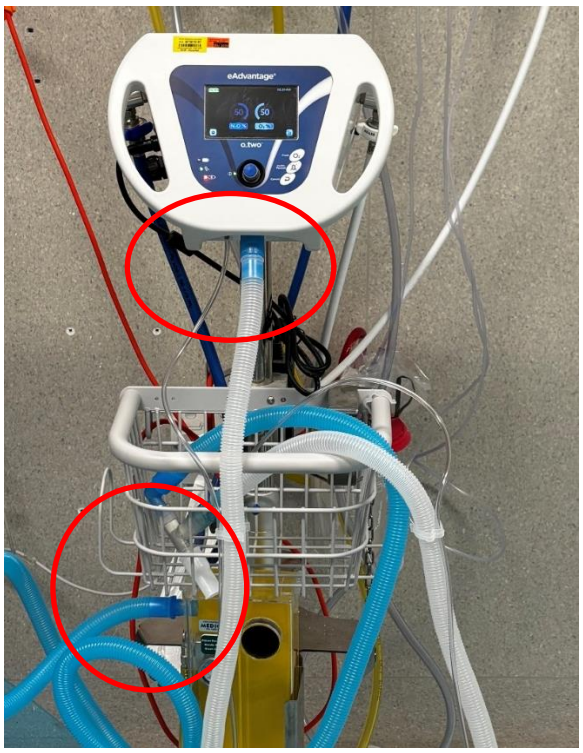
Appendix 1: Demand Administration Setup Using Piped Nitrous Oxide / Oxygen and O-Two eAdvantage® Equipment


For comprehensive diagram of all parts of the Nitrous Oxide O-Two eAdvantage®, refer to Appendix 4.

Ensure the gas hoses (nitrous oxide and oxygen) scavenger tubing and suction regulator are connected to the wall outlets.	
Connect the nitrous oxide breathing circuit with an appropriately sized mask or mouthpiece. Attach the inspiratory tube (white) to the nitrous oxide outlet located on the front of the machine and the expiratory tube (blue) to the scavenging outlet on the scavenger unit.	
Turn on the nitrous oxide demand unit by holding the on/off button.	


Turn on the wall suction used to provide scavenging suction and set the regulator to -15 kPa by ensuring the needle reaches the green indicator on the scavenger unit.	

Appendix 2: Demand Administration Setup Using Portable Nitrous Oxide / Oxygen and O-Two eAdvantage® Equipment

Portable Nitrous Oxide / Oxygen Equipment Setup	
Portable Demand Unit	
Verify that both gas cylinders are adequately filled and connected to the gas hoses. Open the cylinder valves as directed. If a cylinder is empty, replace it before use. There is no required sequence for turning the cylinders on.	

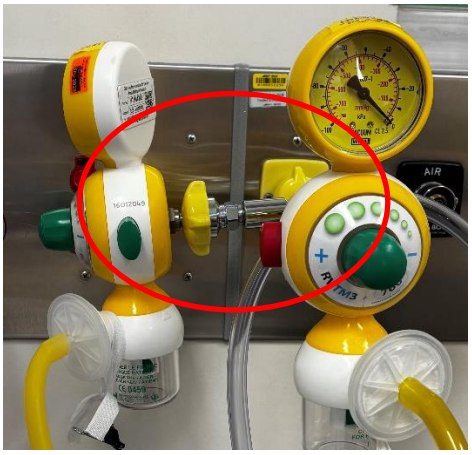Connect the portable suction regulator to the wall suction regulator via the 3-way adapter (circled) to provide scavenging suction.	
Connect the nitrous oxide breathing circuit with an appropriately sized mask or mouthpiece. Attach the inspiratory tube (white) to the nitrous oxide outlet located on the front of the machine and the expiratory tube (blue) to the scavenging outlet on the scavenger unit.	
Turn on machine by holding the on/off button	
Turn on the wall suction used to provide scavenging suction and set the regulator to -15 kPa by ensuring the needle reaches the green indicator on the scavenger unit.	

Appendix 3: Oversedation Treatment

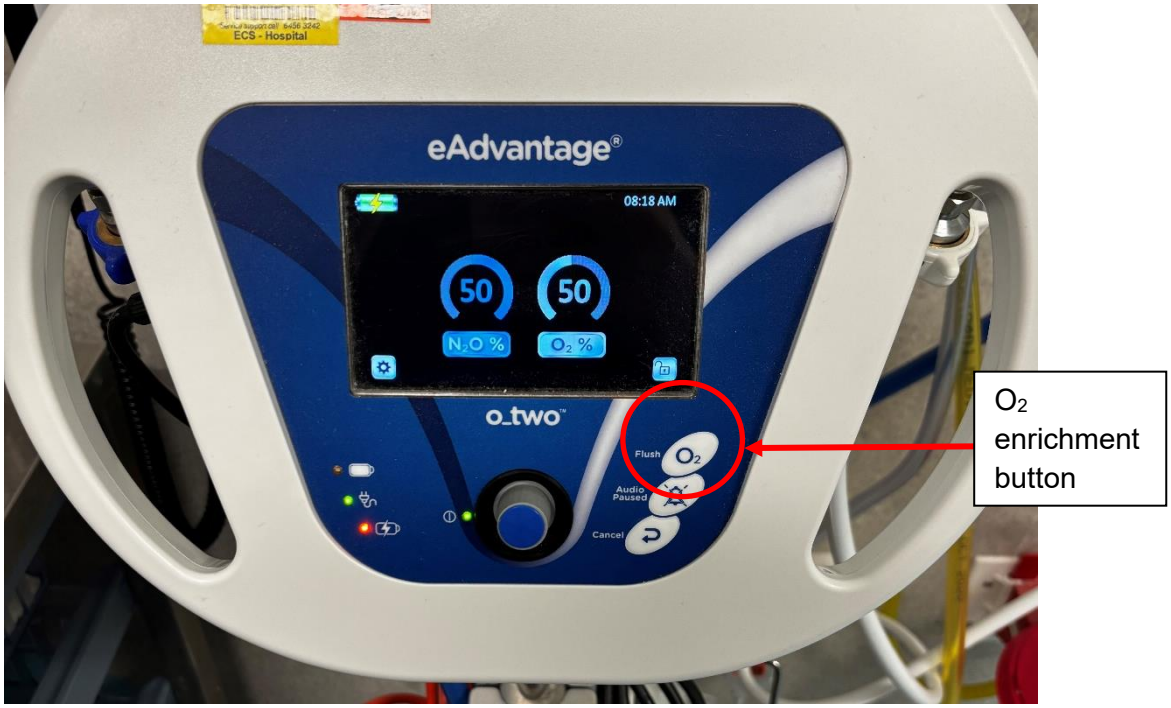Oversedation Treatment
In the event of oversedation, press the 'O_2' enrichment button on the O-Two eAdvantage® as circled in the image below. 

Appendix 4: O-Two eAdvantage® Labelled Diagram

